
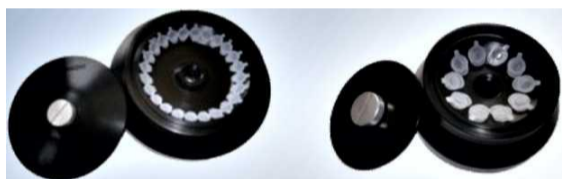


## KT-TGL18E Tabletop High Speed Refrigerated Centrifuge

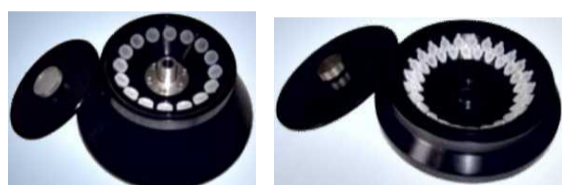

### Features:

1. Adopt carbon-free brush frequency conversion motor, maintenance-free, high torque, fast lifting speed. Fluorine-free compressor unit is adopted, which has no environmental pollution and precise temperature control.
2. Using microcomputer processor for precise control, digital display of speed, temperature, time and other parameters, key programming, switching display operating parameters and RCF value, 10 sets of programs can be stored and recalled, and 10 kinds of lifting rates are available for selection.
3. The independent motor servo of the door cover reduces the operator's labor intensity, automatic protection for over-speed and over-temperature, unbalanced protection, and the body adopts high-quality steel structure.
4. Unique design of spring cone sleeve to connect the rotor and the spindle, loading and unloading the rotor is quick and simple, non-directional, safe and reliable, and easy to use.
5. Multiple rotating heads for selection, can design various adapters according to the test requirements, multi-purpose, three-stage vibration reduction to ensure the optimal centrifugal effect.

### Technical Data

Model	KT-TGL18E/KT-TGL18ER
Max Speed	18500rpm
Max RCF	23810 xg
Max Capacity	6*100ml
Speed Accuracy	± 10rpm
Temperature Range	-20°C ~ +40°C
Temperature Accuracy	± 1°C
Display	LED/LCD
Timer Range	1 -99min/9h59min
Motor	AC brushless motor
Power Supply	AC220V 50HZ 10A
Compressed Units	Imported environmentally friendly refrigerant
Full Power	1100W
Noise	≤62db
Device Dimension	510*630 *380mm(L*W*H)
Net Weight	110kg

**Matched Rotors(customized product available)**

NO	ITEM	MAX SPEED(rpm)	MAX RCF(xg)	VOLUME(ml)	TUBE SIZE(mm)
NO.1	Angle rotor	18500rpm	23810xg	12*1.5/2.2ml	Φ10*42
NO.2	Angle rotor	16000rpm	20653xg	24*1.5/2.2ml	Φ10*42
NO.3	Angle rotor	14000rpm	19450xg	30*1.5/2.2ml	Φ10*42
NO.4	Angle rotor	13000rpm	17930xg	48*1.5/2.2ml	Φ10*42
NO.5	Angle rotor	15000rpm	18510xg	30*0.2/0.5ml	Φ8*30
NO.6	Angle rotor	14000rpm	19970xg	48*0.2/0.5ml	Φ8*33
NO.7	Angle rotor	15000rpm	19350xg	8/10*5ml	Φ14*54
NO.8	Angle rotor	15000rpm	19350xg	16*5ml	Φ14*54
NO.9	Angle rotor	15000rpm	13856xg	3*8*0.2ml (PCR tube)	Φ6*22
NO.10	Angle rotor	14000rpm	12826xg	6*8*0.2ml (PCR tube)	Φ6*22
NO.11	Angle rotor	14000rpm	16370xg	12*7ml	Φ15*62
NO.12	Angle rotor	14000rpm	16460xg	6*10ml	Φ16*87
NO.13	Angle rotor	12000rpm	14856xg	12*10ml	Φ16*87
NO.14	Angle rotor	12000rpm	13215xg	8*15ml	Φ16*118
NO.15	Angle rotor	5000rpm	3515xg	24*15ml	Φ16*118
NO.16	Angle rotor	5000rpm	3830xg	30*15ml	Φ16*118
NO.17	Angle rotor	12000rpm	14510xg	8*20ml	Φ22*104
NO.18	Angle rotor	14000rpm	16050xg	6*30ml	Φ25*94
NO.19	Angle rotor	14000rpm	16050xg	4*50ml	Φ29*106
NO.20	Angle rotor	12000rpm	15635xg	6*50ml	Φ29*106
NO.21	Angle rotor	10000rpm	10813xg	6*70ml	Φ32*150
NO.22	Angle rotor	12000rpm	15635xg	4*100ml	Φ37*124
NO.23	Angle rotor	10000rpm	11380xg	6*100ml	Φ37*124
NO.24	Angle rotor	12000rpm	15869xg	24pcs capillaires	Φ4*100
NO.25	Horizontal ELISA plate rotor	4000rpm	2310xg	2*2*48 hole	-